


## Anti-GOLPH2 Antibody (A10100)

### Specifications:


Name:	Anti-GOLPH2 Antibody
Description:	Rabbit polyclonal antibody to GOLPH2.
Applications:	WB, IHC, ICC/IF
Recommended Dilutions:	WB: 1:500-1:1,000, IHC: 1:50-1:200, ICC/IF: 1:50-1:200
Reactivity:	Human, Mouse, Rat
Immunogen:	Recombinant fusion protein containing a sequence corresponding to amino acids 112-270 of human GOLPH2 (NP_057632.2).
Sequence:	TGERLIRVLQDQLKTLQRNYGRLQQDVLQFQKNQTNLERKFSYDLSQCINQMKEVKEQ CEERIEEVTKKGN EAVASRD LSENNDQRQQLQALSEPQPR LQAAGLPHTTEVPQGKGNV LGNSKSQTPAP SSEVVLDSKRQVEKEETNEIQVVNEEPQRDRL
Host:	Rabbit
Clonality:	Polyclonal
Isotype:	IgG
Conjugate:	Unconjugated
Purification:	Affinity purification.
Molecular Weight:	75 kDa
Product Form:	Liquid
Formulation:	Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol and 0.05% Proclin 300.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.

## Anti-GOLPH2 Antibody (A10100)


### Images:

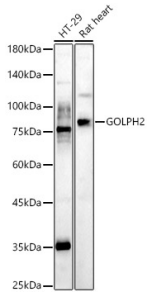


Western blot analysis of various lysates, using Anti-GOLPH2 Antibody (A10100) at 1:1,000 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25 $\mu$ g per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 180s.

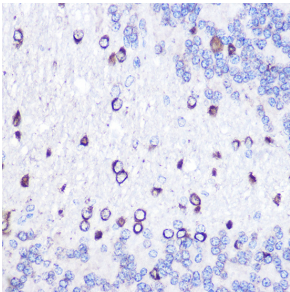

Immunohistochemistry analysis of paraffin-embedded rat brain using Anti-GOLPH2 Antibody (A10100) at a dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

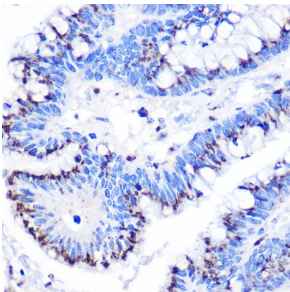

Immunohistochemistry analysis of paraffin-embedded human colon carcinoma tissue using Anti-GOLPH2 Antibody (A10100) at a dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

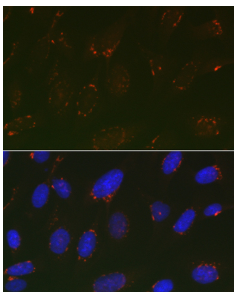

Immunofluorescence analysis of U2OS cells using Anti-GOLPH2 Antibody (A10100) at a dilution of 1:100 (40x lens). DAPI was used to stain the cell nuclei (blue).