


Anti-GALC Antibody (A14482)

Specifications:


Name:	Anti-GALC Antibody
Description:	Rabbit polyclonal antibody to GALC.
Applications:	WB, IHC, ICC/IF
Recommended Dilutions:	WB: 1:500-1:2,000, IHC: 1:50-1:200, ICC/IF: 1:50-1:200
Reactivity:	Human, Mouse, Rat
Immunogen:	Recombinant fusion protein containing a sequence corresponding to amino acids 40-300 of human GALC (NP_000144.2).
Sequence:	GGAYVLDDSDGLGREFDGIGAVSGGGATSRLLVNYPEPYRSQILDYLFKPNFGASLHI LKVEIGGDGQTTDGTPEPSHMHYALDENYFRGYEWWMKEAKKRNPNTLIGLPWSFPG WLGKGFDPYVNLQLTAYYVVTWIVGAKRYHDLDDIDYIGIWNERSYNANYIKILRKML NYQGLQRVKIIASDNLWESISASMLLDAELFKVVDVIGAHYPGTHSAKDAKLTGKKLW SSEDFSTLNSDMGAGCWGRILNQNYINGY
Host:	Rabbit
Clonality:	Polyclonal
Isotype:	IgG
Conjugate:	Unconjugated
Purification:	Affinity purification.
Molecular Weight:	77 kDa
Product Form:	Liquid
Formulation:	Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol and 0.02% Sodium Azide.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.

Anti-GALC Antibody (A14482)


Images:

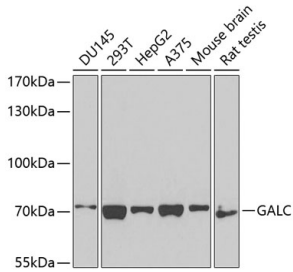


Western blot analysis of extracts of various cell lines, using Anti-GALC Antibody (A14482) at 1:1,000 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25 μ g per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 90s.

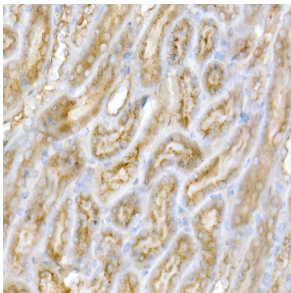

Immunohistochemistry analysis of paraffin-embedded human kidney using Anti-GALC Antibody (A14482) at a dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

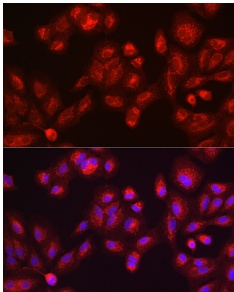

Immunofluorescence analysis of A-549 cells using Anti-GALC Antibody (A14482) at a dilution of 1:200 (40x lens). DAPI was used to stain the cell nuclei (blue).

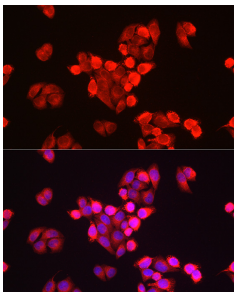

Immunofluorescence analysis of HeLa cells using Anti-GALC Antibody (A14482) at a dilution of 1:200 (40x lens). DAPI was used to stain the cell nuclei (blue).