

## Anti-ZAP70 (phospho Tyr493) Antibody (A16675)

## Specifications:

Name: Anti-ZAP70 (phospho Tyr493) Antibody

Description: Rabbit polyclonal antibody to ZAP70 (phospho Tyr493).

Applications: WB, IHC, ICC/IF

Recommended Dilutions: WB: 1:500-1:1,000, IHC: 1:50-1:100, ICC/IF: 1:50-1:200

Reactivity: Human, Mouse, Rat

Immunogen: A synthetic phosphorylated peptide around Y493 of human ZAP70 (NP\_001070.2).

Sequence: SYYTA

Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Conjugate: Unconjugated

Purification: Affinity purification.

Molecular Weight: 70 kDa

Product Form: Liquid

Formulation: Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol and 0.01% Thiomersal.


Storage: Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.

Disclaimer: This product is for research use only. It is not intended for diagnostic or therapeutic use.




## Anti-ZAP70 (phospho Tyr493) Antibody (A16675)


## Images:

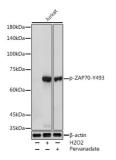

Western blot analysis of extracts of Jurkat cells, using Anti-ZAP70 (phospho Tyr493) Antibody (A16675) at 1:1,000 dilution. Jurkat cells were treated by Hydrogen Peroxide (2 mM) at 37°C for 2 minutes after serum-starvation overnight. Jurkat cells were treated by Pervanadate (1 mM) at 37°C for 30 minutes. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25µg per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 1s.

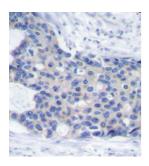

Immunohistochemistry analysis of paraffin-embedded human breast carcinoma tissue using Anti-ZAP70 (phospho Tyr493) Antibody (A16675). Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

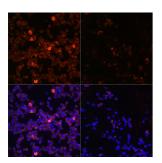

Immunofluorescence analysis of Jurkat cells using Anti-ZAP70 (phospho Tyr493) Antibody (A16675) at a dilution of 1:100 (40x lens). Jurkat cells were treated by Hydrogen Peroxide (2 mM) at 37°C for 2 minutes after serum-starvation overnight. DAPI was used to stain the cell nuclei (blue).