

Anti-Ferritin Light Chain Antibody [FTL/1386] (A248646)

Specifications:

Name: Anti-Ferritin Light Chain Antibody [FTL/1386]

Description: Mouse monoclonal [FTL/1386] antibody to Ferritin Light Chain.

Applications: ELISA, WB, IHC

Recommended Dilutions: WB: 1-2 μg/ml, IHC: 0.25-0.5 μg/ml

Reactivity: Human

Immunogen: Recombinant fragment, around amino acids 38-165, of human FTL protein. The exact

sequence is proprietary.

Host: Mouse

Clonality: Monoclonal

Clone ID: FTL/1386

Isotype: IgG2b

Light Chains: kappa

Conjugate: Unconjugated

Purification: Protein A/G chromatography.

Concentration: 200 µg/ml

Product Form: Liquid

Formulation: Supplied in 10mM Phosphate Buffered Saline with 0.05% BSA and 0.05% Sodium Azide.

Storage: Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.

General Notes: This monoclonal antibody is also available in a different formulation without BSA and

Sodium Azide - Anti-Ferritin Light Chain Antibody [FTL/1386] - BSA and Azide free


(A251828).

Disclaimer: This product is for research use only. It is not intended for diagnostic or therapeutic use.



Anti-Ferritin Light Chain Antibody [FTL/1386] (A248646)


Images:

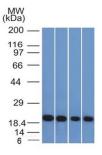

Western blot analysis of A431, HeLa, liver, and testis lysates using Anti-Ferritin Light Chain Antibody [FTL/1386].



Western blot analysis of human HeLa cell lysate using Anti-Ferritin Light Chain Antibody [FTL/1386].

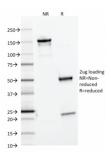

SDS-PAGE analysis of Anti-Ferritin Light Chain Antibody [FTL/1386] under non-reduced and reduced conditions; showing intact IgG and intact heavy and light chains, respectively. SDS-PAGE analysis confirms the integrity and purity of the antibody.



Anti-Ferritin Light Chain Antibody [FTL/1386] (A248646)

Images continued:



Analysis of protein array containing more than 19,000 full-length human proteins using Anti-Ferritin Light Chain Antibody [FTL/1386]. Z-Score and S- Score: The Z-score represents the strength of a signal that a monoclonal antibody (MAb) (in combination with a fluorescently-tagged anti-IgG secondary antibody) produces when binding to a particular protein on the HuProtTM array. Z-scores are described in units of standard deviations (SD's) above the mean value of all signals generated on that array. If targets on HuProtTM are arranged in descending order of the Z-score, the S-score is the difference (also in units of SD's) between the Z-score. S-score therefore represents the relative target specificity of a MAb to its intended target; a MAb is considered to be specific to its intended target, if the MAb has an S-score of at least 2.5. For example, if a MAb binds to protein X with a Z-score of 43 and to protein Y with a Z-score of 14, then the S-score for the binding of that MAb to protein X is equal to 29.