


Anti-CD34 Antibody (A87833)

Specifications:

Name:	Anti-CD34 Antibody
Description:	Rabbit polyclonal antibody to CD34.
Applications:	WB, IHC, ICC/IF
Recommended Dilutions:	WB: 1:500-1:1,000, IHC: 1:50-1:200, ICC/IF: 1:50-1:200
Reactivity:	Human, Mouse, Rat
Immunogen:	A synthetic peptide corresponding to a sequence within amino acids 50-150 of human CD34 (NP_001020280.1).
Sequence:	SNVSTNVSYQETTTPSTLGSTSLHPVSQHGNEATTNITETTVKFTSTSVITSVYGNTN SSVQSQTSVISTVFTTPANVSTPETTLKPSLSPGNVSDLSTTS
Host:	Rabbit
Clonality:	Polyclonal
lsotype:	lgG
Conjugate:	Unconjugated
Purification:	Affinity purification.
Molecular Weight:	105 kDa
Product Form:	Liquid
Formulation:	Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol and 0.01% Thiomersal.
Storage:	Shipped at 4° C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.

Images:

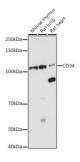


Western blot analysis of extracts of various cell lines, using Anti-CD34 Antibody (A87833) at 1:1,000 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25µg per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Enhanced Kit (RM00021). Exposure time: 180s.


antibodies

Anti-CD34 Antibody (A87833)


Images continued:

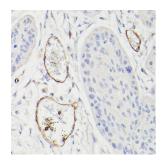

Immunohistochemistry analysis of paraffin-embedded human esophageal cancer using Anti-CD34 Antibody (A87833) at a dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

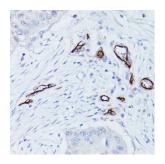

Immunohistochemistry analysis of paraffin-embedded human lung cancer using Anti-CD34 Antibody (A87833) at a dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

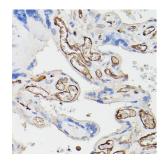

Immunohistochemistry analysis of paraffin-embedded human placenta using Anti-CD34 Antibody (A87833) at a dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.