antibodies

Anti-AMPK alpha 2 Antibody (A92373)


Specifications:

Name:	Anti-AMPK alpha 2 Antibody
Description:	Rabbit polyclonal antibody to AMPK alpha 2.
Applications:	WB, ICC/IF
Recommended Dilutions:	WB: 1:500-1:1,000, ICC/IF: 1:50-1:200
Reactivity:	Human, Mouse, Rat
Immunogen:	Recombinant fusion protein containing a sequence corresponding to amino acids 343-552 of human AMPKa2 (NP_006243.2).
Sequence:	ASSPPSGSFMDDSAMHIPPGLKPHPERMPPLIADSPKARCPLDALNTTKPKSLAVKKA KWHLGIRSQSKPYDIMAEVYRAMKQLDFEWKVVNAYHLRVRRKNPVTGNYVKMSLQLY LVDNRSYLLDFKSIDDEVVEQRSGSSTPQRSCSAAGLHRPRSSFDSTTAESHSLSGSL TGSLTGSTLSSVSPRLGSHTMDFFEMCASLITTLAR
Host:	Rabbit
Host: Clonality:	Rabbit Polyclonal
Clonality:	Polyclonal
Clonality: Isotype:	Polyclonal IgG
Clonality: Isotype: Conjugate:	Polyclonal IgG Unconjugated
Clonality: Isotype: Conjugate: Purification:	Polyclonal IgG Unconjugated Affinity purification.
Clonality: Isotype: Conjugate: Purification: Molecular Weight:	Polyclonal IgG Unconjugated Affinity purification. 62 kDa
Clonality: Isotype: Conjugate: Purification: Molecular Weight: Product Form:	Polyclonal IgG Unconjugated Affinity purification. 62 kDa Liquid


antibodies

Anti-AMPK alpha 2 Antibody (A92373)


Images:

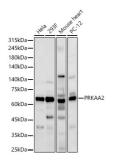


Western blot analysis of various lysates, using Anti-AMPK alpha 2 Antibody (A92373) at 1:900 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25µg per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 10s.

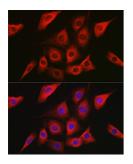

Immunofluorescence analysis of NIH/3T3 cells using Anti-AMPK alpha 2 Antibody (A92373) at a dilution of 1:100 (40x lens). DAPI was used to stain the cell nuclei (blue).

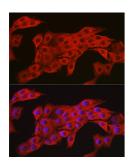

Immunofluorescence analysis of PC-12 cells using Anti-AMPK alpha 2 Antibody (A92373) at a dilution of 1:100 (40x lens). DAPI was used to stain the cell nuclei (blue).

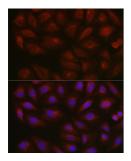

Immunofluorescence analysis of U2OS cells using Anti-AMPK alpha 2 Antibody (A92373) at a dilution of 1:100 (40x lens). DAPI was used to stain the cell nuclei (blue).