


## Anti-STAT3 Antibody (A13365)

### Specifications:

Name:	Anti-STAT3 Antibody
Description:	Rabbit polyclonal antibody to STAT3.
Applications:	WB, IHC, ICC/IF
Recommended Dilutions:	WB: 1:500-1:1,000, IHC: 1:50-1:200, ICC/IF: 1:50-1:200
Reactivity:	Human, Mouse, Rat
Immunogen:	A synthetic peptide corresponding to a sequence within amino acids 600-700 of human STAT3 (NP_003141.2).
Sequence:	TKPPGTFLLRFSSESSKEGGVTFTWVEKDISGKTQIQSVEPYTKQLNNMSFAEIIMGY KIMDATNILVSPLVLYLPDIPKEEAFGKYCRPESQEHPEADPG
Host:	Rabbit
Clonality:	Polyclonal
Isotype:	IgG
Conjugate:	Unconjugated
Purification:	Affinity purification.
Molecular Weight:	86 kDa
Product Form:	Liquid
Formulation:	Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol and 0.05% Proclin 300.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.


### Images:

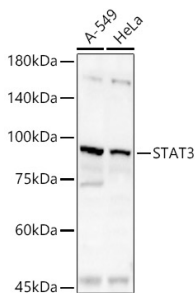


Western blot analysis of various lysates, using Anti-STAT3 Antibody (A13365) at 1:2,500 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25µg per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 60s.

## Anti-STAT3 Antibody (A13365)


Images continued:

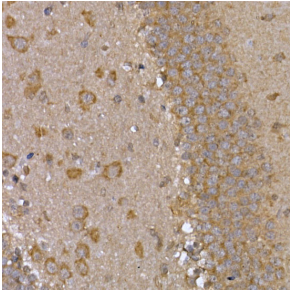

Immunohistochemistry analysis of paraffin-embedded mouse brain using Anti-STAT3 Antibody (A13365) at a dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

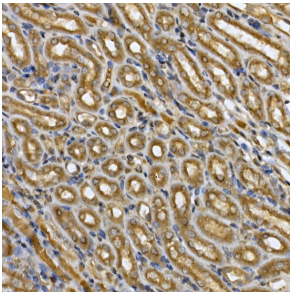

Immunohistochemistry analysis of paraffin-embedded mouse kidney using Anti-STAT3 Antibody (A13365) at a dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

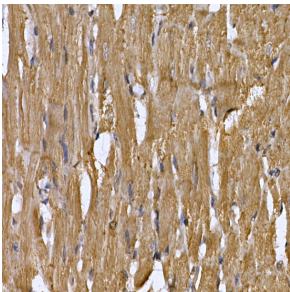

Immunohistochemistry analysis of paraffin-embedded rat heart using Anti-STAT3 Antibody (A13365) at a dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.