

Anti-Aquaporin 3 Antibody (A14218)

Specifications:

Name: Anti-Aquaporin 3 Antibody

Description: Rabbit polyclonal antibody to Aquaporin 3.

Applications: WB, IHC, ICC/IF

Recommended Dilutions: WB: 1:500-1:1,000, IHC: 1:50-1:200, ICC/IF: 1:50-1:200

Reactivity: Human, Mouse, Rat

Immunogen: Recombinant fusion protein containing a sequence corresponding to amino acids 178-292

of human AQP3 (NP_004916.1).

Sequence: IVDPYNNPVPRGLEAFTVGLVVLVIGTSMGFNSGYAVNPARDFGPRLFTALAGWGSAV

FTTGQHWWWVPIVSPLLGSIAGVFVYQLMIGCHLEQPPPSNEEENVKLAHVKHKEQI

Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Conjugate: Unconjugated

Purification: Affinity purification.

Molecular Weight: 32 kDa

Product Form: Liquid

Formulation: Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol and 0.01% Thiomersal.


Storage: Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.

Disclaimer: This product is for research use only. It is not intended for diagnostic or therapeutic use.




Anti-Aquaporin 3 Antibody (A14218)


Images:

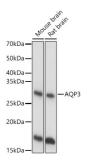


Western blot analysis of extracts of various cell lines, using Anti-Aquaporin 3 Antibody (A14218) at 1:1,000 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25 μ g per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 10s.

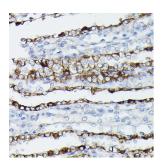

Immunohistochemistry analysis of paraffin-embedded rat kidney using Anti-Aquaporin 3 Antibody (A14218) at a dilution of 1:50 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

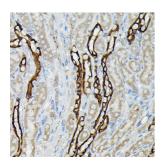

Immunohistochemistry analysis of paraffin-embedded mouse kidney using Anti-Aquaporin 3 Antibody (A14218) at a dilution of 1:50 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

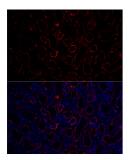

Immunofluorescence analysis of rat kidney cells using Anti-Aquaporin 3 Antibody (A14218) at a dilution of 1:100 (40x lens). DAPI was used to stain the cell nuclei (blue).