

Anti-Somatostatin Receptor 5 Antibody [ARC1282] (A307261)


Specifications:

Anti-Somatostatin Receptor 5 Antibody [ARC1282]
Rabbit monoclonal [ARC1282] antibody to Somatostatin Receptor 5.
WB, IHC
WB: 1:500-1:2,000, IHC: 1:50-1:200
Human, Mouse, Rat
A synthetic peptide corresponding to a sequence within amino acids 1-100 of human SSTR5 (P35346).
MEPLFPASTPSWNASSPGAASGGGDNRTLVGPAPSAGARAVLVPVLYLLVCAAGLGGN TLVIYVVLRFAKMKTVTNIYILNLAVADVLYMLGLPFLATQN
Rabbit
Monoclonal
ARC1282
lgG
Unconjugated
Affinity purification.
45 kDa
Liquid
Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol, 0.05% BSA, and 0.02% Sodium Azide.
Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
This product is for research use only. It is not intended for diagnostic or therapeutic use.


antibodies

Anti-Somatostatin Receptor 5 Antibody [ARC1282] (A307261)


Images:

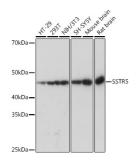

Western blot analysis of extracts of various cell lines, using Anti-Somatostatin Receptor 5 Antibody [ARC1282] (A307261) at 1:1,000 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25µg per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 3s.

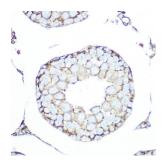

Immunohistochemistry analysis of paraffin-embedded rat testis using Anti-Somatostatin Receptor 5 Antibody [ARC1282] (A307261) at a dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

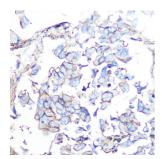

Immunohistochemistry analysis of paraffin-embedded human lung cancer using Anti-Somatostatin Receptor 5 Antibody [ARC1282] (A307261) at a dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.



Immunohistochemistry analysis of paraffin-embedded mouse kidney using Anti-Somatostatin Receptor 5 Antibody [ARC1282] (A307261) at a dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.