

Anti-KIR3DS1 Antibody (A16377)

Specifications:

Name: Anti-KIR3DS1 Antibody

Description: Rabbit polyclonal antibody to KIR3DS1.

Applications: WB, IHC

Recommended Dilutions: WB: 1:500-1:2,000, IHC: 1:50-1:200

Reactivity: Human, Mouse, Rat

Immunogen: Recombinant fusion protein containing a sequence corresponding to amino acids 80-340 of

human KIR3DS1 (NP_001077008.1).

Sequence: FNMSPVTTAHAGNYTCRGSHPHSPTGWSAPSNPMVIMVTGNHRKPSLLAHPGPLVKSG

ERVILQCWSDIMFEHFFLHKEWISKDPSRLVGQIHDGVSKANFSIGSMMRALAGTYRC YGSVTHTPYQLSAPSDPLDIVVTGLYEKPSLSAQPGPKVQAGESVTLSCSSRSSYDMY HLSREGGAHERRLPAVRKVNRTFQADFPLGPATHGGTYRCFGSFRHSPYEWSDPSDPL

LVSVTGNPSSSWPSPTEPSSKSGNLRHLH

Host: Rabbit

Clonality: Polyclonal

Isotype: IgG

Conjugate: Unconjugated

Purification: Affinity purification.

Molecular Weight: 60 kDa

Product Form: Liquid

Formulation: Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol and 0.02% Sodium

Azide.


Storage: Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.

Disclaimer: This product is for research use only. It is not intended for diagnostic or therapeutic use.




Anti-KIR3DS1 Antibody (A16377)


Images:

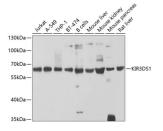

Western blot analysis of extracts of various cell lines, using Anti-KIR3DS1 Antibody (A16377) at 1:1,000 dilution. The secondary antibody was Goat Anti-Rabbit IgG H&L Antibody (HRP) at 1:10,000 dilution. Lysates/proteins were present at 25 μ g per lane. The blocking buffer used was 3% non-fat dry milk in TBST. Detection was with a ECL Basic Kit. Exposure time: 60s.

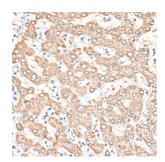

Immunohistochemistry analysis of paraffin-embedded human liver using Anti-KIR3DS1 Antibody (A16377) at a dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

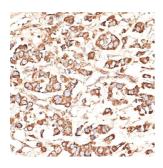

Immunohistochemistry analysis of paraffin-embedded human mammary cancer using Anti-KIR3DS1 Antibody (A16377) at a dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.