


Anti-CD35 Antibody [To5] - BSA and Azide free (A251442)

Specifications:


Name:	Anti-CD35 Antibody [To5] - BSA and Azide free
Description:	Mouse monoclonal [To5] antibody to CD35.
Specificity:	This antibody recognizes a protein of 210-220kDa, which is identified as the complement receptor 1 (CR1)/CD35. This MAb does not block CR1 activity. It is highly specific to CR1 and shows no cross-reaction with CR2. The primary function of CR1 is to serve as the cellular receptor for C3b and C4b, the most important components of the complement system leading to clearance of foreign macromolecules. The Knops blood group system is a system of antigens located on this protein. Follicular dendritic cells (FDC) are restricted to the B-cell regions of secondary lymphoid follicles. They are CD21+/CD35+/CD1a-. This MAb labels follicular dendritic cells and follicular dendritic cell sarcoma.
Applications:	Flow Cytometry, IF, IHC-P
Recommended Dilutions:	Flow Cytometry: 1-2 µg/million cells, IF: 1-2 µg/ml, IHC-P: 1-2 µg/ml
Reactivity:	Human
Immunogen:	Recombinant human CR1 protein.
Host:	Mouse
Clonality:	Monoclonal
Clone ID:	To5
Isotype:	IgG1
Light Chains:	kappa
Conjugate:	Unconjugated
Purification:	Protein A/G chromatography.
Concentration:	1 mg/ml
Product Form:	Liquid
Formulation:	Supplied in 10mM Phosphate Buffered Saline; without Sodium Azide and carrier free.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
General Notes:	This monoclonal antibody is also available in a different formulation with BSA and Sodium Azide - Anti-CD35 Antibody [To5] (A248260).
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.

Anti-CD35 Antibody [To5] - BSA and Azide free (A251442)

Images:

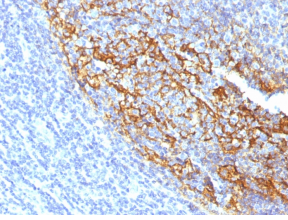

Immunohistochemical analysis of formalin-fixed, paraffin-embedded human tonsil using Anti-CD35 Antibody [To5].

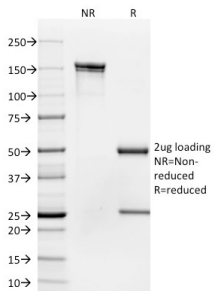

SDS-PAGE analysis of Anti-CD35 Antibody [To5] under non-reduced and reduced conditions; showing intact IgG and intact heavy and light chains, respectively. SDS-PAGE analysis confirms the integrity and purity of the antibody.