


Anti-ZAP70 Antibody [ZAP70/2047] - BSA and Azide free (A253526)

Specifications:


Name:	Anti-ZAP70 Antibody [ZAP70/2047] - BSA and Azide free
Description:	Mouse monoclonal [ZAP70/2047] antibody to ZAP70.
Specificity:	ZAP70 is a 70kDa protein tyrosine kinase found in T-cells and natural killer cells. Control of this protein translation is via the IgVH gene. ZAP70 protein is expressed in leukemic cells of approximately 25% of chronic lymphocytic leukemia (CLL) cases as well. Anti-ZAP70 expression is an excellent surrogate marker for the distinction between the Ig-mutated (anti-ZAP70 negative) and Ig-unmutated (anti-ZAP70 positive) CLL subtypes and can identify patient groups with divergent clinical courses. The anti-ZAP70 positive Ig-unmutated CLL cases have been shown to have a poorer prognosis.
Applications:	ELISA, Flow Cytometry, IF, IHC-P
Recommended Dilutions:	Flow Cytometry: 1-2 µg/million cells, IF: 1-2 µg/ml, IHC-P: 1-2 µg/ml
Reactivity:	Human
Immunogen:	Recombinant fragment, around amino acids 247-382, of human ZAP70 protein. The exact sequence is proprietary.
Host:	Mouse
Clonality:	Monoclonal
Clone ID:	ZAP70/2047
Isotype:	IgG1
Light Chains:	kappa
Conjugate:	Unconjugated
Purification:	Protein A/G chromatography.
Concentration:	1 mg/ml
Product Form:	Liquid
Formulation:	Supplied in 10mM Phosphate Buffered Saline; without Sodium Azide and carrier free.
Storage:	Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
General Notes:	This monoclonal antibody is also available in a different formulation with BSA and Sodium Azide - Anti-ZAP70 Antibody [ZAP70/2047] (A250346).
Disclaimer:	This product is for research use only. It is not intended for diagnostic or therapeutic use.

Anti-ZAP70 Antibody [ZAP70/2047] - BSA and Azide free (A253526)


Images:

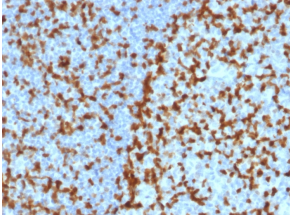


Immunohistochemical analysis of formalin-fixed, paraffin-embedded human tonsil using Anti-ZAP70 Antibody [ZAP70/2047].

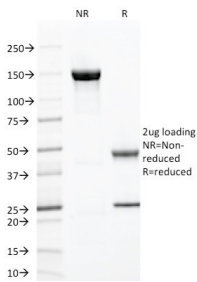

SDS-PAGE analysis of Anti-ZAP70 Antibody [ZAP70/2047] under non-reduced and reduced conditions; showing intact IgG and intact heavy and light chains, respectively. SDS-PAGE analysis confirms the integrity and purity of the antibody.

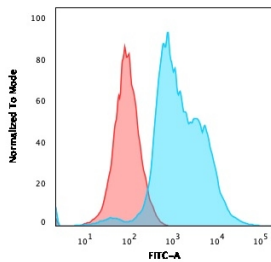

Flow cytometric analysis of PFA fixed Jurkat cells using Anti-ZAP70 Antibody [ZAP70/2047] followed by Goat Anti-Mouse IgG (CF® 488) (Blue). Isotype Control (Red).

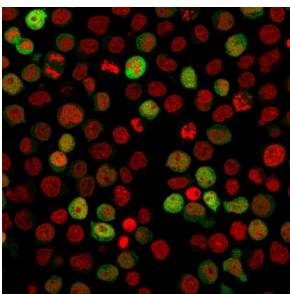

Immunofluorescent analysis of paraformaldehyde fixed Jurkat cells stained with Anti-ZAP70 Antibody [ZAP70/2047] followed by Goat Anti-Mouse IgG (CF® 488) (Green). Nuclear counterstain is RedDot.