

# Anti-MDM2 (phospho Ser166) Antibody [ARC2909] (A309403)

#### Specifications:

Name: Anti-MDM2 (phospho Ser166) Antibody [ARC2909]

Description: Rabbit monoclonal [ARC2909] antibody to MDM2 (phospho Ser166).

Applications: IHC

Recommended Dilutions: IHC: 1:50-1:200

Reactivity: Human, Mouse, Rat

Immunogen: A phospho synthetic peptide corresponding to residues surrounding S166 of Human MDM2.

Host: Rabbit

Clonality: Monoclonal

Clone ID: ARC2909

Isotype: IgG

Conjugate: Unconjugated

Purification: Affinity purification.

Product Form: Liquid


Formulation: Supplied in Phosphate Buffered Saline, pH 7.3, with 50% Glycerol, 0.05% BSA, and 0.02%

Sodium Azide.

Storage: Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.

Disclaimer: This product is for research use only. It is not intended for diagnostic or therapeutic use.

#### Images:

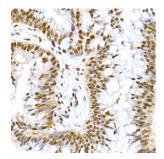


Immunohistochemistry analysis of paraffin-embedded human colon carcinoma tissue using Anti-MDM2 (phospho Ser166) Antibody [ARC2909] (A309403) at a dilution of 1:20 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.



## Anti-MDM2 (phospho Ser166) Antibody [ARC2909] (A309403)


### Images continued:

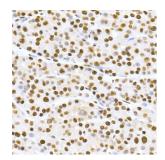

Immunohistochemistry analysis of paraffin-embedded human liver cancer using Anti-MDM2 (phospho Ser166) Antibody [ARC2909] (A309403) at a dilution of 1:20 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.



Immunohistochemistry analysis of paraffin-embedded mouse liver using Anti-MDM2 (phospho Ser166) Antibody [ARC2909] (A309403) at a dilution of 1:20 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

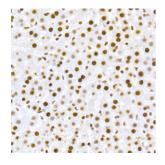

Immunohistochemistry analysis of paraffin-embedded rat liver using Anti-MDM2 (phospho Ser166) Antibody [ARC2909] (A309403) at a dilution of 1:20 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.